

Viking™ M mobile lift

The all-round lift



The Viking M mobile lift is ideal for meeting the more challenging patient handling demands, yet remains easy to manoeuvre, even in confined spaces. Despite its low weight, this durable aluminum lift has a lifting capacity of 205 kg.

Features and benefits

Informative control box displays **number of lift cycles, battery level, lift overload and service indicator**

Optional **quick-release hook system** which means caregivers can easily and rapidly change a sling bar or install a scale

Optional armrests help **support patients during ambulation** as well as making it easier for caregivers to manoeuvre the lift

Three variable height settings accommodate a variety of different lifting situations


Intuitive and ergonomically-designed hand control incorporates indicators for low battery level and when the Safe Working Load of the lift is exceeded

An extensive assortment of accessories is available including slings, slingbars and scales

Viking M mobile lift



Hand control

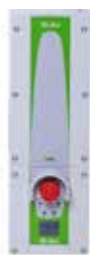

Built-in battery and control box

Technical Specifications

Maximum load	205 kg
Weight	33 kg
Height	1340 - 2020 mm
Protection class	IPX4
Hand control protection class	IPX4
Slingbar options	
standard	Universal 450 Slingbar 450 mm
optional	Mini Slingbar 220 Universal Slingbar 350 Universal TwinBar 670 Crossbars
Base	
height	105 mm
length	1230 mm
width	690 - 1110 mm
expansion	electric
Castors	
front	75 mm double band
rear	75 mm double band
Material	aluminium
Emergency lifting	electric
Emergency lowering	manual and electric

Lifting speed (without load)	36 mm / sec
Lifting interval	1270 mm
Maximum lifting height	1790 mm
Sound level	46 dbA
Battery	lead-acid gel
Charger	built-in
Control Box	
audible low battery signal	when <25% power left
service information	displayed
lift counter	displayed
overload indicator	displayed*
battery power indicator	displayed*
Hand control	wired
Product number	2040045

*also displayed on the hand control

Hill-Rom is a leading global medical technology company with approximately 10,000 employees worldwide. We partner with health care providers in more than 100 countries by focusing on patient care solutions that improve clinical and economic outcomes in five core areas: **Advancing Mobility, Wound Care and Prevention, Clinical Workflow, Surgical Safety and Efficiency, and Respiratory Health.** Around the world the people, programs, and products of Hill-Rom work towards one mission: **Enhancing outcomes for patients and their caregivers.**

Not all products/options are available in all countries. For further information about our products and services, please contact your local Hill-Rom representative or visit our webpage:

www.hill-rom.com

Viking M is a versatile lift intended mainly for use in health care, intensive care and rehabilitation. Class I
Manufacturer: Liko AB, Nedre vägen 100, SE-975 92 Luleå, Sweden

This medical device is a regulated health product which, pursuant to such regulation bears a CE mark. Hill-Rom recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2016 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.
Doc. No: 5EN183307-02, 01 August 2016

Enhancing outcomes for patients and their caregivers:

Hill-Rom